

## Safety

*\*Refer to the manufacturer's manual or website for more information about your equipment.*

- + Use a seatbelt if needed to maintain proper position.
- + Keep weight centered on the chair when reaching for objects.
- + Keep anti-tippers in down position.
- + Do not use wheelchair to do "wheelies".
- + Do not use wheelchair on an escalator.
- + Do not use the wheelchair for sitting in a vehicle.
- + Use the wheel locks to stabilize chair.
- + Keep hands and fingers away from moving parts of chair.

## How to Use

### Unfolding

1. Place your palm on the seat rail and press downward until the chair is fully open. Be sure your hand placement doesn't allow for your fingers to get caught.
2. Engage both wheel locks and open the legrests.


### Folding

1. Swing the legrests to the front, locked position.
2. Push the footplates to the vertical position.
3. Grasp of middle of the seat at the front and back edges and pull up to close the frame.

### Transferring

1. Remove or swing away armrest and/or legrest that may be in the way, if possible.
2. Position chair as close as possible to seat where you are transferring.
3. Lock the wheels.
4. Shift body weight to the seat.

*\*Only try transfers by yourself if you have good mobility and upper body strength. Practice transferring with a caregiver. Use a transfer board if needed.*

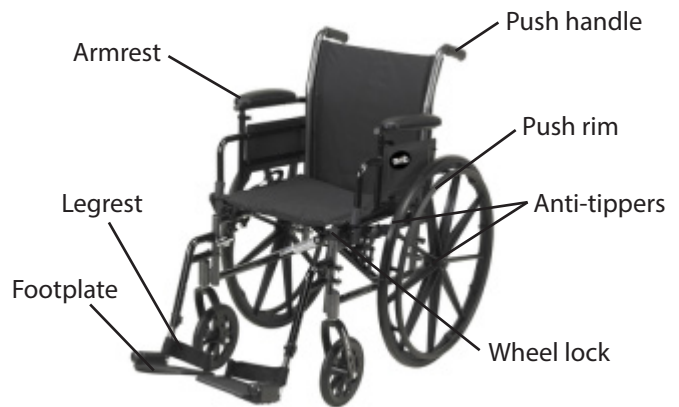

## For Caregivers

### Going forward over a curb

1. Tilt the chair back to balance point and roll the chair forward until the front casters clear the curb.
2. Set the casters down and push the chair forward so the rear wheels and anti-tippers roll up and over the curb.

### Going backward over a curb

1. Turn the anti-tippers to the up position.
2. Stand on the upper sidewalk with the rear wheels against the curb.
3. Adjust the anti-tippers back to down position.
4. Pull the chair up and over the curb.

## Equipment Care

### Cleaning

- + Wipe down metal parts and vinyl upholstery. Use a soft cloth and a mild cleaner.
- + Use upholstery cleaner to clean upholstery when needed.

### Wear & Tear

- + Examine tires and caster, and look for loose rims and cracks.
- + Check the spokes for proper position and connection to wheel.
- + Look for wobbling motion in wheels.